

A RANGE OF PRODUCTS THAT  
FACILITATE MORE EFFECTIVE SURGERY



A commitment to the highest quality and an ability to exceed international standards is what makes Frankenman and our Chex range of products stand out. We put health first. This manifests itself in our total dedication to a zero defect outcome. Our manufacturing, production, research and development, and customer service facilities are the envy of the international medical device sector.

We source and use only the highest quality materials and components from around the world. This means that not only do we adhere to internationally recognized measures of quality, such as CE and ISO 13485 certification, we regularly exceed

of our staff is second to none and our design capabilities are of the highest standards throughout the industry culminating in the award of full ISO 9001: 2000 certification for our Quality Management systems.




WHEN YOU PUT PATIENTS' HEALTH FIRST,  
REMARKABLE THINGS FOLLOW



State-of-the-art manufacturing facilities.

INFORMED DESIGN.  
SCRUPULOUSLY APPLIED.  
METICULOUSLY PRODUCED.

PROFESSOR ULRICH BAUMGARTNER OF THE KREISKLINIK  
OTTOBUREN PRESENTS HIS EXPERIENCE OF THE USE OF THE  
CHEX CPH RANGE OF STAPLERS

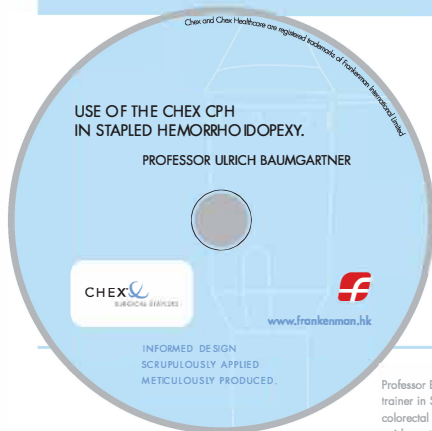

PROFESSOR ULRICH  
BAUMGARTNER

CHIEF OF THE DEPARTMENT OF  
GENERAL, VISCERAL SURGERY AND  
COLO-PROCTOLOGY SPECIALIST FOR  
SURGERY AND VISCERAL SURGERY,  
PROCTOLOGY (SPECIALIST AREAS OF  
INTEREST: COLORECTAL SURGERY)

Professor Baumgartner is a recognised and experienced  
trainer in Stapled Hemorrhoidopexy, STARR and other  
colorectal procedures and has produced an excellent video  
guide on the techniques he employs.

### OUR RANGE OF PRODUCTS INCLUDE

SINGLE USE CIRCULAR  
STAPLER FOR RECTAL  
PROLAPSE AND  
HAEMORRHOIDS



SINGLE USE  
RELOADABLE LINEAR  
STAPLER



SINGLE USE CURVED  
INTRALUMINAL  
CIRCULAR STAPLER



CPH34HV



CS COMPACT™



SINGLE USE  
TRANSVERSE CUTTING  
LINEAR STAPLER



CIRCULAR STAPLERS FOR RECTAL  
PROLAPSE AND HEMORRHOIDS

INFORMED DESIGN.  
SCRUPULOUSLY APPLIED.  
METICULOUSLY PRODUCED.

## A RANGE OF PRODUCTS TO MEET THE NEEDS OF MODERN SURGICAL PRACTICE

The introduction of the Chex™ CPH 34 and CPH 32 range of Circular Staplers for Stapled Hemorrhoidopexy has been received with the highest praise from Surgeons around the world. Undoubtedly...

"suma cum laude"

## A RANGE OF PRODUCTS THAT, QUITE SIMPLY, ARE NOW CONSIDERED THE BEST

When you put health first, you test everything as strenuously as possible. You test it yourself and have others test it for you, you ask questions, you go back if necessary, you incorporate what you've learned into the manufacturing process.

You devote all your energies to product quality and end-user cost reduction. Only then can you be satisfied that the products you produce are truly the best possible.

At Frankenman, we believe we're not so much in the business of making instruments but in the business of making surgeries more successful.

That's why we listen to doctors who've put our instruments to use. We read their reports. We invite them to think with us, to brainstorm, to help us develop new, innovative products.

That way we know we will be able to deliver the highest quality to the medical profession both today and way into the future.

### CPH 34 ACTUAL SIZE

21.60cm<sup>3</sup> INTERNAL HOUSING VOLUME

TOP VIEW

1

SMOOTH, ERGONOMIC "UP FREE" CONICAL DESIGN OF ANVIL

2

LONGER "ANVIL ROD" OF ONE PIECE CONSTRUCTION ENABLING BETTER ACCESS AND VISUALIZATION WHEN SECURING "PURSE STRING" AND GREATER MECHANICAL STABILITY FACILITATING ACCURATE STAPLE FORMATION

3

FIRST TO INTRODUCE A MINIMUM TISSUE GAP OF 0.75MM AS STANDARD

4

CPH 34 FEATURES 4 SUTURE CHANNELS PROVIDING MORE OPTIONS FOR ACCURATE AND TENSION-FREE PLACEMENT OF PULSURES

10

"WING NUT" TO ASSIST CONTROL AND DEXTERITY WHEN OPENING AND CLOSING ANVIL

### CPH 32 ACTUAL SIZE

17.17cm<sup>3</sup> INTERNAL HOUSING VOLUME

SIDE VIEW

1

2

3

4

5

6

7


8

ONLY PRODUCT THAT HAS 32 STAPLES IN 2 ROWS TO ASSIST MORE EFFECTIVE HEMOSTASIS

HIGHEST QUALITY SURGICAL STEEL CUTTING BLADE

UNIQUE "CUTTING" MECHANISM WITH CIRCUMFERENTIAL PUSH ASSEMBLY

RADICAL REDESIGN OF STAPLE HOUSING GIVING SURGEON ACCURATE CONTROL OF TISSUE DRAWN INTO HOUSING  
LOWERED TISSUE STOP TO FACILITATE SMOOTH, CONTROLLABLE ENTRY OF TISSUE INTO SUTURE HOUSING

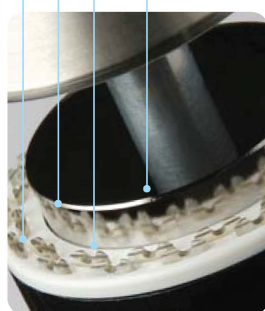

### ACCESSORIES

9

FIRST TO SUPPLY TWO ANAL DIATORS IN ADDITION TO STANDARD ANAL DIATOR. ALL PRODUCTS ARE SUPPLIED WITH "BUTTERFLY" DIATOR ASSISTING WHEN ACCESS IS RESTRICTED



### STAPLE INFORMATION

MODEL	CARTRIDGE COLOUR	STAPLE ROWS	OUTSIDE DIAMETER	CUTTING DIAMETER	STAPLE LEG LENGTH	STAPLE CROWN LENGTH	STAPLE COUNT	PER BOX
CPH 34	GREEN	2	34.5MM	26MM	4.2MM	3.8MM	32	3
CPH 32	WHITE	2	32.5MM	24MM	4.2MM	3.8MM	32	3

### KEY FEATURES

- LARGE INTERNAL HOUSING CAPACITY TO ACCOMMODATE MORE TISSUE
- ADDITIONAL STAPLES TO PROVIDE EFFECTIVE HEMOSTASIS
- WING NUT CLOSING HANDLE FACILITATES EASE OF CLOSURE

FDA CE 0123

**REY MEDICAL SA: YOUR PARTNER IN SURGERY**

• QUALITY • SECURITY • ECONOMY

Visit our website: [www.reymedical.ch](http://www.reymedical.ch) or call us: 022 743 24 10