

FRANKENMAN

CS STANDARD PORTFOLIO





WHEN YOU PUT PATIENTS' HEALTH FIRST, REMARKABLE THINGS HAPPEN



FRANKENMAN are very proud of our Circular stapler portfolio. Simply the most extensive of any manufacturer on the market.

Many devices have been designed with the input of Global Key opinion leaders to bring unique iterations and innovative devices to the surgical community whilst providing improved functionality, ergonomics and the continual evolution of our circular stapler portfolio, embodies our ethos of "When you put patients' health first, remarkable things happen"



SAFETY & EFFICACY

Longer, 5.0mm, open staple leg height, to ensure correct "B" Shaped formation & security,

even in thicker tissue.

Unique engineering of integrated staple driver & 4 cutting blade, gives greater stability & strength leading to improved performance during the 4 firing cycle.

Robust Safety locking mechanism prevents premature firing.

Clear indicators for opening & closing.

Extended pursestring tying area, reduces risk of inaccurate application.

Reduced "force to fire" vs alternative devices. due to design of washer in anvil, reduces distal movement, & increases stability during the firing cycle. This may also reduce the risk of incomplete firing.




ANVIL HOLDER

- All Frankenman Circular Staplers are supplied with an anvil holding rod as standard.
- Allows secure handling of the anvil when applying a pursestring suture.
- **C** Controlled handover and docking of the anvil during TransAnal procedures.



VARIABLE TISSUE COMPRESSION

- **f** Ensures appropriate staple formation for a variety of tissue thicknesses.
- "White", "Blue" & "Green" staple closure heights in one device.
- Gives the surgeon greater control.

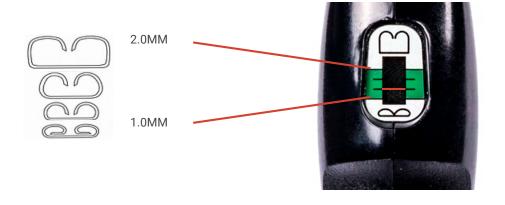

PRODUCT CODE	DESCRIPTION		TOTAL WORKING LENGTH	OUTER DIAMETER (mm)	LUMEN SIZE (mm)	STAPLE QUANTITY	OPEN STAPLE HEIGHT	CLOSED STAPLE HEIGHT	STAPLE CROWN	PER BOX
CS21L	LONG ENDOSTEER™		520MM	21	13	16	4.5MM	1.0 - 2.0MM	1 4MM	3
CS25L*	LONG ENDOSTEER™			25	17	20	5.0MM			
CS28L*	LONG ENDOSTEER™			28	20	24				
CS32L	LONG ENDOSTEER™			32	24	28				
CS21	STANDARD		420MM	21	13	16	4.5MM	2.0		
CS25	STANDARD			25	17	20	5.0MM			
CS28*	STANDARD			28	20	24				
CS32*	STANDARD			32	24	28				

^{*}References available in our stock at REY Médical SA




QUALITY CONTROL

We pride ourselves on world leading continual quality improvement processes.

ALL stapling devices and reloads have a unique serial number, as well a batch/lot number, to allow for more accurate tracking and traceability.

In addition we test fire 100% of our products (before they are then reloaded and packaged), to verify correct B-Shaped closure.

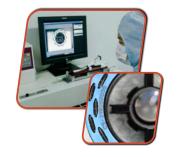

YOUR PARTNER IN SURGERY SINCE 2003

INNOVATION AND QUALITY AT THE SERVICE OF SURGERY

+41 (0)22 743 24 10 • contact@reymedical.ch • www.reymedical.ch

The best technology available on the market, at the best price, at the right time